

Peer Grading Rubric PHSX 516 Effects of var. HV and emission

	Reviewer Author	Kelly Title of p	aper Ma
Ounline of I	Good	Fair	Poor
Quality of da	to improve results.	ed Completed all aspect of experiment but no attempt to cross che	Followed manual,
Title	Evokes interest in physics	Plain statement of what was done	Misleading
Abstract	Clearly states what was done, provides final result & explains significance	States what was don but with no physics context	e Does not include fi result nor put the result into context
Introduction	Sets current experiment in physics context of the time. Explains theory or theories that are relevant N <sub>citations</sub> > 2	Explains theory but does not highlight significance of experiment N <sub>citations</sub> =	Poor explanation of theory N <sub>citations</sub> = 0
Method  Tust use Fig 2 of Fig 7	Reader feels that they really understand how to repeat the experiment. Diagram of apparatus is clear.	Most steps are explained but some missing. Figure may be poor.	Explanation is garble Figure is missing or very poor.
Results  data columns into another to be pietled?	This section begins with "Figure 2 shows stopping voltage "Figures are clear and uncluttered. Legend font at least 12pt Each fig cited. Systematic & statistical uncertainties estimated.	Figure are cited but for some Figures no conclusions are made about the data. Captions are incomplete. Legends smaller than 12pt text.	text. Captions are unclear, Figures are
Incertainties ummary	X <sup>2</sup> is reasonable Recapitulates what was measured and why it matters.	X² unreasonable States final result but does not put it in	No uncertainty analysi Missing or makes claims not supported
riting	Compares to other results Spelling & grammar correct. Paper gradually builds up an argument-based on data.	A few grammar mistakes. Text doesn't clarify what can be earned from data.	by data.  Many grammar mistakes. Text does not make a coherent argument.
	inal report addresses all	Some corrections	Final report essentially the same as draft

Something that world be interested to see plotted would Zmin VS voltage, I think it would be linear

be copied format

٧



E=hv

E = hc = V = hc  $= \sqrt{\lambda}$   $h = \sqrt{\lambda}$ 

# Effects of variable high voltage and emission current on X-ray energy spectrum of molybdenum.

Thresa Kelly and Michael Murray
University of Kansas, PHSX 616 Physical Measurements
(Dated: October 5, 2023)

## ABSTRACT

#### I. INTRODUCTION

The atomic modal has gone through several major changes throughout history. John Dalton proposed the "solid sphere" model of the atom in his 1803. He claimed that atoms are indivisible particles characterized by their mass that can be combined in whole-number ratios to make compounds. Next, in 1897, Joseph John Thomson discovered a negatively-charged particle—the electron. Thus, the atom was not indivisible; instead, it had a complex internal structure. This supported Lord Kelvin's model, which claimed that the atom was a sphere of positive charge with pockets of negatively-charged electrons. The atomic model was further refined in 1911 by Ernest Rutherford. From his famous gold-foil experiment, Rutherford concluded that atoms must have a small, dense, positively charged core-the nucleus-with the electrons orbit randomly around it. Soon after, Niels Bohr proposed his quantized shell model. Here, the electrons orbit the nucleus of fixed energy. Electrons can absorb energy from a photon and jump to a higher orbital with a higher energy level, or emit energy and jump to a lower level.

Bremsstrahlung radiation occurs when highly energetic, free electrons decelerate, moving to a lower energy level, when they contact matter. This produces a continuous spectrum of X-ray photons. Electrons penetrating deep into the atomic shell produce characteristic emission lines that peak above the continuous spectrum. These eject the innermost orbital electrons via collisions, which create gaps that the outer-orbital electrons fall into. Energetic, free electrons can be produced by an electrical voltage or current. Electrons emitted from the cathode are incident on the anode, and produce X-ray emission lines unique to the anode material.

One method of observing an X-ray spectra is by reflecting the X-rays off of a crystal. Bragg's law of refraction (Equation 1) relates the scattering angle  $(\theta)$  to the wavelength  $(\lambda)$  of first order diffraction [1]. In this research, we use a NaCl crystal with a lattice plane spacing (d) of 282.01 pm.

$$\lambda = 2d\sin\theta \tag{1}$$

The energy of the X-ray radiation can be calculated using the following electromagnetic frequency equations: Equation 2 (where  $\nu$  is the frequency and c is the speed of light)

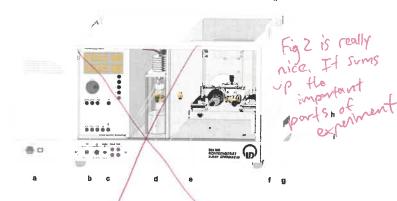


FIG. 1. 554 800 Rontgengerat X-ray apparatus, which consists of the following components: (a) main power panel, (b) control panel, (c) connection panel, (d) X-ray tube chamber, (e) experiment chamber with goniometer, (f) fluorescent screen, (g) free, (h) lock, (i) feet, and (k) carrying handles [3].

$$\nu = \frac{c}{\lambda} \tag{2}$$

and Equation 3 (where h is Planck's constant) [2].

$$E = h\nu \tag{3}$$

In this research, we investigate the x-ray energy spectrum of molybdenum X-ray tube. We observe the effects of varying the high voltage and emission current on the resulting spectrum.

#### II. METHOD

For this research, we use the 554 800 Rontgengerat X-ray apparatus (see Figure 1). The apparatus is controlled by a microprocessor, which controls the X-ray tube, goniometer, and Geiger-Muller counter tube [3]. The X-ray tube is a hot cathode tube with a molybdenum (Mo) anode. The anode material is  $K_{\alpha}=17.4$  keV and  $K_{\beta}=19.6$  keV with a voltage range of 0.0 kV < U < 35 kV and current range of 0.0 mA < I < 1 mA. The goniometer consists of two independent stepper motors to control the position of the target stage and sensor arm. It has an angular resolution of 1° across a range of  $-10^{\circ}$  to 170°. The NaCl crystal is placed on the target stage. The target and sensor arm can be coupled to maintain

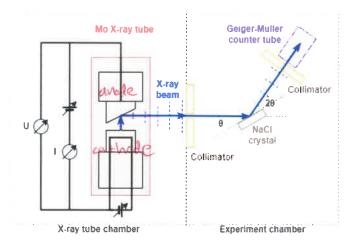


FIG. 2. Experiment setup. The X-ray beam (shown in blue) is generated by the X-ray tube with molybdenum (Mo) anode (red). The diffracted X-ray beam is collimated (orange) as it enters the experiment chamber. The beam reflects off the NaCl (gray) at an angle  $\theta$  before passing through a second collimator (orange) and being absorbed by the Geiger-Muller counter tube. The counter tube is coupled to be at an angle of  $2\theta$  with reference to the horizontal X-ray beam. The emission current (I) and high voltage (V) are variable.

a 2:1 angle. The Geiger-Muller counter tube is a sensor that can count the number of incident X-rays, up to 9999 s<sup>-1</sup>. The control panel allows the user to set the X-ray tube high voltage (U), emission current (I), measuring time per angular step  $(\delta t)$ , goniometer angular step width  $(\delta \beta)$ , and goniometer upper and lower limits  $(\beta)$  [3]. The NaCl crystal in Bragg configuration in combination with the counter tube comprise a spectrometer.

In this research, we are investigating the effects of the high voltage and emission current on the X-ray energy spectrum by first order Bragg reflection on an NaCl crystal. To do this, we conduct two experiments: (a) varying the high voltage at a constant current, and (b) varying the emission current for a constant voltage. The setup for these experiments is shown in Figure 2. For both experiments, the slit opening of the X-ray tube collimator is set to be 5 cm from the center of the NaCl crystal on the target stage; the Geiger-Muller sensor is then adjusted to be 6 cm from the crystal. The measuring time is set to  $\delta t = 10$  s, angular step width  $\delta \beta = 0.1^{\circ}$ , and goniometer upper limit of  $\beta = 2.5^{\circ}$  and lower limit  $\beta = 12.5^{\circ}$ . The apparatus is also set to  $2\theta$  coupling between the target and sensor on the goniometer. For experiment (a), we take measurements between the goniometer limits for each U = 15, 20, 25, 30, and 35 keV while the maintaining a constant I = 1.00 mA. For experiment (b), we vary take measurements for each  $I=0.40,\,0.60,\,0.80,\,$  and 1.00 mA with a constant U=35 kV. We perform two trials for each experiment to cross-check the data and improve results.

	Voltage [±0.1 kV]
$U_0$	15.0
$U_1$	20.0
$U_2$	25.0
$U_3$	30.0
$U_4$	35.0

TABLE I. Five different voltages for the constant I=1.00 mA. The X-ray photon counts are tracked from  $(2.5\pm0.1)^{\circ} \le \beta \le (12.5\pm0.1)^{\circ}$  at  $(0.1\pm0.1)^{\circ}$  resolution.

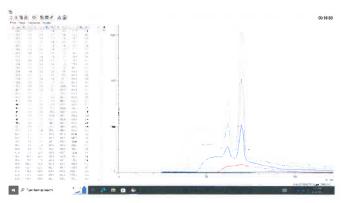


FIG. 3. Trial 1 for constant current and variable voltage. The largest curves show the larger voltages.

#### III. RESULTS AND DISCUSSION

In this research, we investigate the effects of varying voltage and current on the X-ray spectrum for a molybdenum (Mo) anode. We claim the following instrumental uncertainties: I is  $\pm 0.01$  mA, U is  $\pm 0.1$ ,  $\beta$  is  $\pm 0.1$ , and t is  $\pm 1$ .

Table I shows the five voltages used for the constant current of I=1.00 mA. We complete two independent trials, measuring the X-ray spectrum for each U twice. The results are shown in Figures 3 and 3. We observe very good agreement between the two trials. The figures show that increasing the high voltage causes the Bremsstrahlung radiation spectra and emission lines to become more pronounced. The emission features are difficult to resolve at  $U<25~\rm kV$ , suggesting that the X-rays are unable to penetrate the inner orbitals.

### IV. CONCLUSION

<sup>[1]</sup> W. H. S. Bragg and W. L. S. Bragg, The reflection of x-rays by crystals, Nature 91, 477 (1913).

<sup>[2]</sup> M. Planck, Ueber das gesetz der energieverteilung im normalspectrum, Annalen der Physik 309, 553–563 (1901).

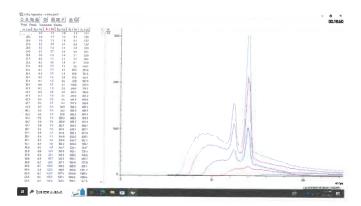


FIG. 4. Trial 2 for constant current and variable voltage. The largest curves show the larger voltages.

[3] L. D. GmbH, Instruction sheet 554 800 (2014).

= 5 · 2		